

DIGITAL RADIOGRAPHY



VENUS IRAY TETHERED DR PLATE

Venu 1717XV is an innovative a-Si Portable Detector for veterinary radiographic imaging. It features a 139 μm pixel pitch, direct deposition CsI scintillator, and 16 bit AID conversion to provided excellent image quality. Venu 1717XV supports a fast work flow for animal X-ray imaging, and is the optimal choice for both retrofit and new veterinary DR system solutions.


Features

- **Cassette detector per ISO 4090, fits in bucky**
- **Direct deposition CsI for better image quality at lower dose**
- **139 μm pixel pitch, with 16 bit ADC for more image details**
- **Multiple trigger modes**
- **Highly reliable, flexible tetherSupports a fast and efficient workflow**

[Read More](#)

CODE E-GEN1131-1

Categories: [Digital Radiography](#), [X-Ray Equipment](#)

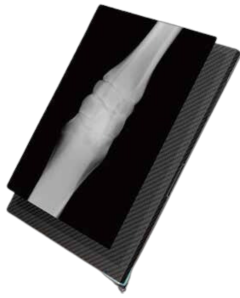

NON GLASS DIGITAL FLAT PANEL DETECTOR

First Flexible TFT CSI X-ray Detector Specifically Designed for Equine Imaging! High-Performance Flexible Wireless a-Si Flat Panel Detector

[Read More](#)

CODE E-GEN1120

Categories: [Digital Radiography](#), [Digital Radiography](#), [Equine](#), [X-Ray Equipment](#)

LONGTAIL- THE FULL SPINE FLAT PANEL DR PACKAGE



The Longtail is the world's first monolithic full-spine flat-panel DR detector. No more long scan wait, no image stitching, no image gaps, and no plate overlap.

It is the first of its kind to produce long-bone and full-spine digital studies instantaneously with the highest geometric accuracy.

At 44 inch (112 cm) long, the Longtail can take a full spine 42 x 17 inch (108 x 43 cm) X-ray images in just one single shot. With tripled power of image acquisition, Longtail reads out images as fast as any regular 14 x 17 inch or 17 x 17 inch detectors and processes outstanding images with high DQE and high resolution. Equipped with highly sensitive AED function, the Longtail can be easily connected to and synchronized with any kind of generator.

Use of the Longtail delivers the lowest dosage for the most sensitive patients, or those who require multiple images over time. It can also reduce dosage for full-spine and long-bone studies by at least 50% while providing the highest geometric accuracy available. The monolithic panel reduces unevenness and provides a perfect density gradient with no overlap.

Longtail DR features:

- The Longtail is fitted with two rechargeable lithium ion batteries that take 2-3 hours to fully charge. One full charge yields 650 exposures and lasts 11 hours in standby mode. For wired use, the Longtail comes with a 5m magnetically connected cable.
- The Longtail detector can also be used for localized studies by simply collimating to the area of clinical interest. The Longtail app will auto-shutter to generate images of any part of the body.
- DROC is a DR image acquisition workstation based on Windows. It provides complete control of all image capture functions for DR application. It delivers the better image quality with lower dose. DROC also has fully integrated with most of imaging components as X-ray generators, flat panel detectors (FPD), collimators, DAPs, and mechanical positioning systems. Intuitive User Interface and Workflow.
- Advanced Clinical applications.

Lead time: Subject to stock availability, transit times and radiation licensing. Estimated shipping date to be provided once the purchase is confirmed

[Read More](#)

CODE E-GEN1098

Category: [Digital Radiography](#)

Tags: [Digital radiography](#), [DR](#), [Full Spine](#), [Imaging Stitching](#), [Leg length](#), [Longtail](#), [Orthopedic Imaging](#), [Spine](#)



ATX ULTRA 99V DR PACKAGE

Highest resolution, 99 micron, accurate imaging.

The ULTRA-99V is a cassette-sized, next generation flat panel DR detector. Battery stands by for at least 12 hours and supports IEEE 802.11ac with transmission speed up to 1300Mbps. The detector acquires outstanding images with high DQE and high resolutions. Equipped with highly sensitive AED function, the detector can be easily connected to and synchronized with any kind of X-ray generator

Included on package:

- ATX Flat DR Panel - 14x17inch - Wireless CsI scintillator
- Lenovo ThinkCentre M920z All-in-One Computer - Core i5 i5-8500 - 8 GB RAM - 512 GB SSD - 60.5 cm (23.8") 1920 x 1080 - Desktop
- VET DROC Acquisition Software
- Training & Installation
- iTX Cloud Pacs DICOM Image Storage platform- one year free trial
- 2 Years standard warranty

Lead time: Subject to stock availability, transit times and radiation licensing. Estimated shipping date to be provided once the purchase is confirmed

[Read More](#)

CODE E-GEN1102

Category: [Digital Radiography](#)

Tags: [14x17](#), [17x17](#), [ATX](#), [Digital radiography](#), [DR Panel](#)





OWANDY OPTEO DR SENSORS (USB)

Designed to stay with you for the long term, the Opteo digital x-ray sensor is sturdy and resistant. The cable connecting it to the USB interface unit has been specially reinforced to offer greater resistance to impacts and twisting.

Easy to use, Opteo offers the best image quality/price ratio. High-quality images with minimal “noise” make for a fast, exact diagnosis every time.

Size 1: 25,2mm x 36,9mm

Size 2: 30,8mm x 43,2mm

Sensors can be supplied with a computer at extra cost. Please contact ATX on **1800 330 118** if you would this option.

[Read More](#)

CODE N/A

Categories: [Digital Radiography](#), [X-Ray Equipment](#)

Tags: [Dental Sensor](#), [Dental xray](#), [Opteo](#), [Opteo DR](#), [X-ray](#)

ATX ULTRA-140V DR PACKAGE

The ULTRA-140V is a cassette-sized, next generation flat panel DR detector. Battery stands by for 8 hours and supports IEEE 802.11ac with transmission speed up to 1300Mbps. The detector acquires outstanding images with high DQE and high resolutions. Equipped with highly sensitive AED function, the detector can be easily connected to and synchronized with any kind of X-ray generator

Included on Package:

- ATX Flat DR Panel - 10x12inch, 14x17inch, 17x17inch - Wireless or Wired / Csl or Gadox
- Lenovo ThinkCentre M920z All-in-One Computer - Core i5 i5-8500 - 8 GB RAM - 512 GB SSD - 60.5 cm (23.8") 1920 x 1080 - Desktop
- Ecom DROC Acquisition Software
- Training & Installation
- ATX Cloud PACS DICOM Image Storage platform
- 2 Years standard warranty

[Read More](#)

CODE E-GEN1101

Categories: [Digital Radiography](#), [X-Ray Equipment](#)

Tags: [14x17](#), [17x17](#), [ATX](#), [Digital radiography](#), [DR Panel](#)



NMI FLAT PANEL DR PACKAGE

Saturn 8000 17x17 is NMI's new slim portable flat panel detector for digital radiography with wide active area of 17" x 17". It features NMI's unique wireless communication method, AP mode allows users directly connect to the detector with laptops, tablet PCs and even smart phones. Furthermore, reliable automatic X-ray exposure detection enables integration with any kind of X-ray systems. Both Csl and GADOX scintillator types are available, with both wired and wireless configurations for each scintillator type.* Wide active area of 17" x 17" * Slim and portable * High spatial resolution with 140um pixel array* Wi-Fi data transfer with dual band (2.4GHz and 5GHz)* Stable and reliable automatic exposure detection* Direct communication with smart devices* Viewer software running on Windows™ OS (DxWorks)

Saturn DR Features:

- Speedy and high-quality image acquisition program
- Acquisition of images exposed from Saturn 8000 Series
- Third-party integration Systems: Generator, U-arm, Collimator, DAP
- DICOM 3.0 compatible: Modality Worklist (MWL), Modality Performed Procedure Step (MPPS), DICOM Send / Print, Storage Commitment, Multi-task patient
- Image acquisition and management for multiple patients at the same time are available
- PureImpact™: post-processing algorithm for an efficient diagnosis with finer details

Lead time: Subject to stock availability, transit times and radiation licensing. Estimated shipping date to be provided once the purchase is confirmed

[Read More](#)

CODE E-GEN1099

Category: [Digital Radiography](#)

Tags: [17x17](#), [DR Panel](#), [NMI](#)





KONICA MINOLTA CR - REGUS SIGMA II

COMPLETE SYSTEM WITH ALL-IN-ONE TOUCH SCREEN COMPUTER AND IMAGE PILOT SOFTWARE.

The REGIUS Σ is Konica Minolta's smallest and lightest CR. The desktop CR has a footprint of only 0.31m² and weighs a mere 28kg.

- **Environmentally Friendly:** The REGIUS Σ was designed with the environment in mind. Energy consumption has been reduced to 100VA - about the same as a house hold light bulb. That's one tenth the power of the smallest film processor. Important with power costs soaring.
- **High Resolution:** The REGIUS Σ supports pixel sizes of 87.5 μ m and 175 μ m for all cassette sizes, achieving a maximum resolution of 4020x4892 pixels (35x43cm at 87.5 μ m).
- **Integral Processing:** Integral Processing - the result of Konica Minolta's vast imaging experience - is a robust algorithm that produces the optimum diagnostic image automatically. This patented technique greatly simplifies x-ray image acquisition while delivering incredibly consistent image quality and virtually eliminating retakes.

[Please contact us for pricing](#)

[Read More](#)

CODE 2KMRegius CR

Category: [Digital Radiography](#)

